



národní  
úložiště  
šedé  
literatury

## **Ride (on Varian) or How to Improve Baseline in $^{29}\text{Si}$ NMR Spectra of Water Glasses**

Schraml, Jan  
2013

Dostupný z <http://www.nusl.cz/ntk/nusl-153271>

Dílo je chráněno podle autorského zákona č. 121/2000 Sb.

Tento dokument byl stažen z Národního úložiště šedé literatury (NUŠL).

Datum stažení: 24.05.2024


Další dokumenty můžete najít prostřednictvím vyhledávacího rozhraní [nusl.cz](http://nusl.cz) .

## RIDE (ON VARIAN) OR HOW TO IMPROVE BASELINE IN $^{29}\text{Si}$ NMR SPECTRA OF WATER GLASSES

Jan Schraml<sup>\*a</sup>, Peter Sandor<sup>b</sup>, Stefan Korec<sup>c</sup> and Bronislav Foller<sup>d</sup>

<sup>a</sup> Institute of Chemical Process Fundamentals of the ASCR, v. v. i., Rozvojová 135, 165 02 Prague 6, Czech Republic. <sup>b</sup> Agilent Technologies R&D a. Marketing GmbH & Co. KG, 76337 Waldbronn, Germany. <sup>c</sup> Vodní sklo a.s., U Tonasa 172 / 2, 403 31 Ústí nad Labem, Czech Republic. <sup>d</sup> University of West Bohemia, New Technologies - Research Centre in the West Bohemian Region, Univerzitní 8, 306 14 Pilsen, Czech Republic.

It is shown by experiments that replacing one-pulse sequence by RIDE (Ring Down Elimination) [1,2] pulse sequence dramatically improves the baseline of  $^{29}\text{Si}$  NMR spectra and eliminates the signal from some probes.

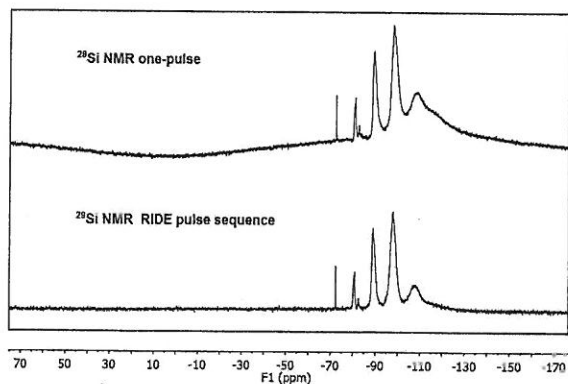

*Acknowledgments:* The authors appreciate the help of Dr. Radek Pohl, Institute of Organic Chemistry and Biochemistry AS CR v. v. i. and Prof. Radovan Fiala of CEITEC. Collaboration with CEITEC within the Open Access Project ID number LM2011020 is gratefully acknowledged. This work was supported by the project FR-TI1/335 "Geopolymeric composites with high technical parameters" provided by the Ministry of Industry and Business of the Czech Republic.

### References

- 1 P. S. Belton, I. J. Cox, R. K. Harris, *J. Chem. Soc. Faraday Trans. 2* **1985**; *81*, 63-75.
- 2 W. Kozmiński, K. Jackowski, *Magn. Reson. Chem.* **2000**; *38*, 459-463.